

Human CA-1600 4 vials 6 chromosomes each
Mouse CA-1601 3 vials 7 chromosomes each

Introduction

The FISH protocol is divided into two stages. Denaturation and Hybridisation are performed on **Day One**. Washing and Detection are performed on **Day Two**.

On day one, the DNA of the chromosomes and paints are denatured and the hybridisation process (reannealing) takes place overnight. On day two, the slides are washed to remove unbound DNA sequences followed by detection, counterstaining and mounting.

Kit Contents		
Product Code	Description	Volume
CA-1600/1	Probe: 4 vials human (Chromosomes 1-6, 7-12, 13-18 and 19-Y) 3 vials mouse (Chromosomes 1-7, 8-14 and 15-Y)	5 tests
1124-F1-50	Detect F1	20µl
1124-F2-50	Detect F2	35µl
1124-DT-50	Detergent (Tween 20)	1ml
1124-Y2-50	Detect Y2 (Cy5 Streptavidin)	20µl

Requirements (not provided)

Equipment	Reagents
Ethanol cleaned slides	Sodium Chloride
Coverslips	Sodium Citrate
Eppendorf tubes	HCl
Coplin jars	Pepsin
Humidified chamber	Formamide
Micro-pipette 1µl, 10µl, 500µl	Absolute Ethanol
Pipette 10ml, 20ml	Deionised distilled water
Vortex	Fixogum rubber cement
Parafilm	Clear nail varnish
Micro-centrifuge	DAPI II®
65°C, 45°C Water bath	
37°C Incubator	
Fluorescence microscope with a suitable filter set	

Approx time:

Slide preparation 90 min
 Denaturation and hybridisation: 30 min + 1 day

This product is for research use only

Cambio Ltd, The Irwin Centre, Scotland Road, Dry Drayton, Cambridge, UK, CB23 8AR
 October 2009

Denaturation and Hybridisation – Day One

Solutions to be prepared: 20XSSC
2XSSC
Denaturation solution
Pepsin solution (optional)

Solution 20XSSC: 87.6g NaCl
44.1g Na Citrate
up to 500ml Deionised Distilled water

Adjust pH to 7.0 using concentrated HCl (before finalising water volume), aliquot and autoclave. Store at 4°C).

Solution 2XSSC: 50ml 20XSSC
450ml Deionised distilled water
500ml 2XSSC.

Denaturation solution: 70ml Formamide
30ml 2XSSC
100ml Denaturation solution.

Pepsin solution: 500µl Stock pepsin solution (1% in water)
49.5ml 10mM HCl
50ml Pepsin solution

Pepsin stock solution can be stored at -20°C in small aliquots.

Note: *Ensure all solutions are mixed well.*

Procedure:

Slide preparation and denaturation:

1. Prepare 4ml of fresh 3:1 fixative (methanol:glacial acetic acid).
2. Place tube of metaphase solution on ice.
3. Use pre-cleaned slides stored in 100% ethanol, wipe off ethanol on the surface with tissue paper and air-dry. Place on a horizontal surface.
4. Resuspend the cells in the fixative by flicking tube gently. Place a single drop of cell suspension on a slide. Check the spreading of chromosomes under phase contrast microscope. If under-spread, add a drop of fresh fixative immediately after cell suspension is placed on slide. Leave slide to dry at room temperature.

Note: *If cytoplasm is still present on the slides, treat the slides with pepsin (see optional pepsin pre-treatment).*

5. Dehydrate by serial ethanol washing for 2 min each in 70% (v/v) ethanol, 70%, 90%, 90%, and 5 min 100%. Dry at room temperature.

This product is for research use only

Cambio Ltd, The Irwin Centre, Scotland Road, Dry Drayton, Cambridge, UK, CB23 8AR

October 2009

Pepsin pre-treatment (optional):

Slides can be pepsin treated on the previous day.

Dehydrate slide in 100% ethanol for 5 min and dry at room temperature.

Incubate slide in pepsin solution for 2-5 min.

(The incubation time will vary depending on amount of cytoplasm present).

Wash in 50ml 2XSSC for 5 min. Repeat and then rinse briefly in distilled water.

Dehydrate by serial ethanol washing for 2 min each in 70% (v/v) ethanol, 70%, 90%, 90%, and 5 min in 100%. Dry at 65°C for 1 hour and continue with Step 7 of the protocol.

6. Age for 60 min on 65°C hot plate. At this stage prepare and denature probe.

Note: *If you have pre-treated the slides a day in advance and aged them overnight at room temperature you do not need to age for 60 min on 65°C hot plate.*

7. Denature slide by incubating in pre-warmed Denaturation solution at 65°C for 1½ -2 min.
8. Quench slides in ice-cold 70% (v/v) ethanol for 4 min and dehydrate by serial ethanol washing for 2 min each in 70% (v/v) ethanol, 70%, 90%, 90%, and 5 min 100%. Dry at room temperature.

Note: *Denaturation of the slides is an important step. Be sure that Denaturation solution is at the right temperature. Some slides benefit from 1½ min denaturation and others up to 2 min. The right timing, which is determined by trial and error, depends on type of cells used, metaphase preparation, brand of formamide, etc.*

Probe preparation and denaturation:

9. Warm probe to 37°C, vortex and centrifuge for 1-3 seconds. Aliquot 10µl probe into a 500µl eppendorf tube.
10. Denature at 65°C-72°C for 10 min and then allow the probes to pre-anneal at 37°C for 10-60 min. The denaturation of probe can either be performed in a water bath, or, ideally, on a PCR machine.

Hybridisation:

11. Apply probe (10µl) onto the slide. Apply coverslip and remove air bubbles by gently pushing on coverslip with a pencil. Seal with rubber cement.
12. Place slide in an air tight, pre-warmed humidified chamber and incubate in the dark at 37°C overnight.

This product is for research use only

Cambio Ltd, The Irwin Centre, Scotland Road, Dry Drayton, Cambridge, UK, CB23 8AR

October 2009

Washing and Detection – Day Two

Solutions to be prepared:

1XSSC
4XSSC
Detergent wash solution
Stringency wash solution
Diluted Detection reagents

Solution 1XSSC: 25ml 20XSSC
 475ml Deionised distilled water
 500ml 1XSSC.

Solution 4XSSC: 100ml 20XSSC
 400ml Deionised distilled water
 500ml 4XSSC.

Detergent wash solution: 500ml 4XSSC
 250µl Detergent (Tween 20)
 500ml Detergent wash solution.

Stringency wash solution: 50ml Formamide
 50ml 1XSSC
 100ml Stringency wash solution.

Stringency wash solution can be reused 5 times but should be discarded after 2 months.

Note: *Ensure all solutions are well mixed.*

Diluted Detection Reagents (DDR):

DDR Y2 (Cy5) + F1 (FITC) 1.20µl Y2
 3.00µl F1
 620.8µl Detergent Wash Solution
 625µl Diluted Detection Reagent (Y2 & F1)

DDR F2 (FITC) 3.0µl F2
 622.0µl Detergent Wash Solution
 625µl Diluted Detection Reagent (F2)

Note: *All volumes of these working solutions are for five slides.*

This product is for research use only

RainbowFISH™ Protocol



Procedure:

Washing:

1. Pre-warm to 45°C in a water bath at least 30 min before starting:
 - Two Coplin jars of Stringency wash solution (50ml each)
 - Three Coplin jars of Solution 1XSSC (50ml each)
 - One Coplin jar of Detergent wash solution (50ml)

Note: *The temperature is important. Check the temperature of the solutions in the Coplin jar and not of the water in the water bath.*

2. Take out the slide from the incubator and carefully remove rubber cement. Place in solution 1XSSC to remove the coverslip.
3. Stringency washes:
 - Wash slides twice by incubating 5 min each in Stringency wash solution (45°C).
 - Wash slides twice by incubating 5 min each in Solution 1XSSC (45°C).
 - Incubate slide for 4 min in Detergent wash solution. (45°C).

Apply 125µl of pre-mixed Detection Reagent Y2 and F1 onto the slide and cover with Parafilm immediately.












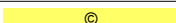




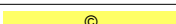


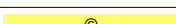



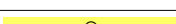

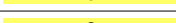






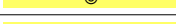
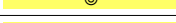

4. Incubate slide in a humidified chamber for 15-20min at 37°C.
5. Remove Parafilm from the slide and wash 3 times 4 min each in Detergent wash solution at room temperature by emptying and refilling the Coplin jar.
6. Apply 125µl Diluted Detection Reagent F2. Cover with parafilm immediately and incubate at 37°C for 15-20 min.
7. Remove Parafilm from the slide and wash 3 times 4 min each in Detergent wash solution at room temperature by emptying and refilling the Coplin jar.
8. Drain slide well and mount with 50µl of DAPI II®.
9. Apply glass coverslip and seal with nail varnish. Store slides in the dark at 4°C.
10. View slides using epifluorescence filters specific for Cy3, Cy5, FITC and DAPI II®.



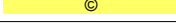



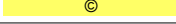























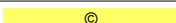


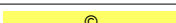

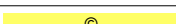
Note: *Capture pictures in the following order: Cy3, Cy5, FITC and DAPI II® and visualise.*

This product is for research use only

Cambio Ltd, The Irwin Centre, Scotland Road, Dry Drayton, Cambridge, UK, CB23 8AR

October 2009

Human CA-1600				
	Human Chromosome	FITC	Cy3	Biotin
Pool I	1			
	2			
	3			
	4			
	5			
	6			
Pool II	7			
	8			
	9			
	10			
	11			
	12			
Pool III	13			
	14			
	15			
	16			
	17			
	18			
Pool IV	19			
	20			
	21			
	22			
	X			
	Y			

Mouse CA-1601				
	Mouse Chromosome	FITC	Cy3	Biotin
Pool I	1			
	2			
	3			
	4			
	5			
	6			
	7			
Pool II	8			
	9			
	10			
	11			
	12			
	13			
	14			
Pool III	15			
	16			
	17			
	18			
	19			
	X			
	Y			

This product is for research use only

Cambio Ltd, The Irwin Centre, Scotland Road, Dry Drayton, Cambridge, UK, CB23 8AR

October 2009